

For more
information,
visit

www.cynosure.com/picosure-pro



¹ RealSelf Worth It Rating as of May 2021.

©2021 Cynosure, LLC. All rights reserved. Cynosure, LLC owns exclusive rights to photography. Use of photography without written permission of Cynosure is prohibited. Cynosure, PicoSure and the Platinum Focus lens are registered trademarks and trademarks of Cynosure, LLC.

These before and after photographs are of a patient treated with Cynosure's PicoSure device, not the PicoSure Pro device advertised herein. The PicoSure Pro device has the same parameters as the PicoSure device and so while individual results vary, similar results are expected.

Like all medical procedures, not all patients are suitable for the treatment. A qualified practitioner is solely responsible for evaluating each subject's suitability to undergo treatment and for informing those being treated about any risks involved with the treatment, pre-and postoperative care, and any other relevant information. Individual results may vary and are not guaranteed.

921-7051-010

Specs & Accessories

Nominal Wavelength	755nm
Pulse Width	500-900 ps
Pulse Energy	300 mJ ± 20%
Max Fluence	6.37 J/cm ²
Spot Sizes	Zoom Handpiece: 2.0 mm-6.0 mm Fixed Handpieces: 5.0mm, 6.0 mm, 8.0 mm, 10.0 mm
Power Supply	200-240 V~, 4.5 kVA, 50/60 Hz, Single Phase
Electrical Service Requirements	30 Amp Dedicated Outlet
Maximum Thermal Emissions	15,500 Btu/h
Dimensions (H x W x D)	43 in (109 cm) x 22 in (56 cm) x 42 in (107 cm)
Weight	375 lb (171 kg)

PicoSure Pro 755 with the Flat Lens is FDA cleared to treat pigmented lesions in skin types I-VI. PicoSure Pro 755 with the Platinum Focus Lens is FDA cleared to treat acne scars and wrinkles in skin types I-IV. Patient results will vary.

PICOSURE[®] *pro*
by CYNOSURE[™]



The art of
PERFORMANCE



Introducing *PicoSure Pro* /

Removes unwanted pigment in *all skin types*

Reduces wrinkles, fine lines, acne scars and pore size with *no downtime*

Improves overall skin tone & texture; providing patients with a lighter, *brighter and more radiant complexion*

The *only* picosecond device FDA-cleared to treat melasma pigment, nevus of Ota and Hori's nevus

10 minute treatment/

Platinum Focus Lens with **40% more shots**

4 fixed handpieces with **adjustable fluence** (*new 5mm handpiece*)

1 Zoom handpiece | 2-6mm

Optional **532nm** and **1064nm** wavelengths for tattoo removal capabilities

New Turbo mode for hard-to-treat tattoos

755 nm.
10 minutes.
Unrivaled results.


Skin revitalization with *no downtime*



Unleash the artist in you with *PicoSure Pro* /

You offer a variety of treatments, each carefully researched and selected with your patients' best interests in mind. Passion for the art and science that encompasses skin revitalization treatments gives you an artist's eye when it comes to providing your patients with the results they desire.

Our unrivaled scientific research and engineering of the newest picosecond laser technology in the market will enable you to take the art of skin revitalization to new levels while providing your patients' with the best outcomes possible.

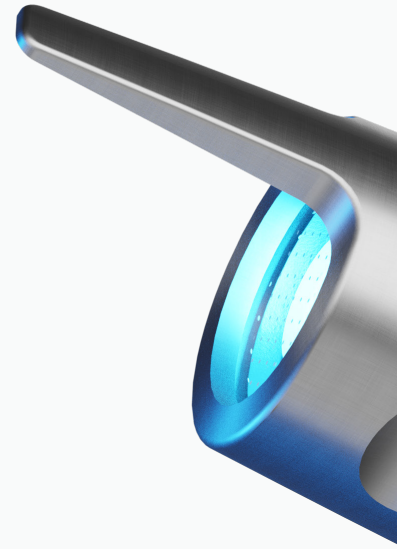



Utilizes 755nm energy to optimally treat skin revitalization with no downtime

40% more shots so that you can treat more patients

New adjustable fluence capability for added customization

Revolutionizing the industry in picoseconds/

First to Market

Cynosure invented the world's first 755nm picosecond laser, revolutionizing the aesthetics industry

Safe & Clinically Proven

More clinical evidence than any other picosecond device
>83 peer reviewed clinical publications

Best in Class Reliability & Service

Manufactured in the U.S. with the highest quality parts, backed by world-class engineering and best-in-class practice support

Globally Trusted Brand

Over a million treatments successfully performed

History of *Innovation* ✓

Cynosure's quest to develop an efficacious skin revitalization treatment with less pain and downtime is backed by extensive clinical research and innovation.

1



Acolade
2008

2



Affinity
2009

3



Medlite & Revlite
2011

4



PicoSure
2012



5



PicoSure Pro | 2022
The newest advancement
in picosecond technology

Platinum *Focus*[™] Lens


Treats fine lines,
wrinkles, acne scars
and large pores with
no downtime

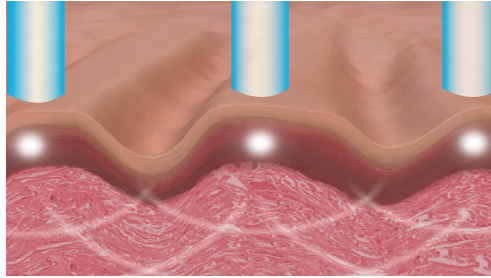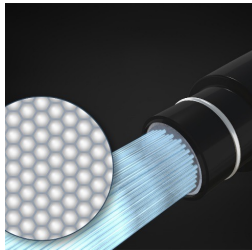

Low intensity and high intensity energy create an elegant injury in the epidermis called intraepidermal Laser Induced Optical Break-downs or LIOBs. As the LIOBs is formed, it creates a pressure wave that expands into the dermis.



Histology shows the formation of LIOBs in the epidermis.

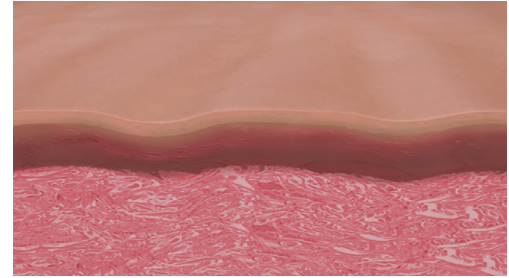



The pressure wave activates the body's natural healing response to increase collagen and elastin production without the thermal side effects or downtime of other lasers.




Histology shows the dermal response which results in the production of new collagen and elastin.

Flat Lens

Treats unwanted pigmentation in all skin types

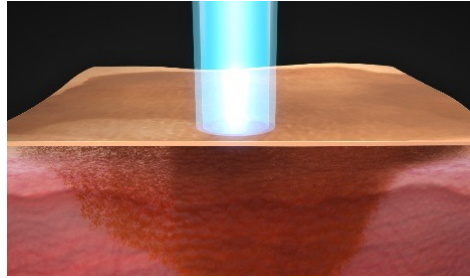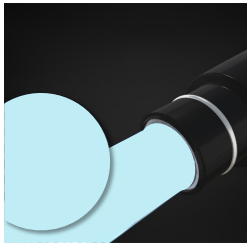

Effectively shatters unwanted pigment into tiny sand-like particles which are naturally absorbed and removed by the body.

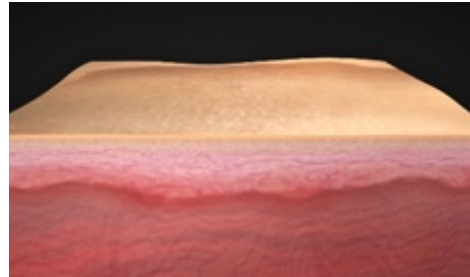

Visibly clears epidermal and dermal pigment without harming the surrounding skin.






The Power of 755nm/

Treat all your patients'
skin needs *safely* &
effectively

Effective skin revitalization with no downtime

The 755 wavelength has 3x the absorption ratio of melanin to blood (vs. 1064). This makes it the ideal wavelength suited to revitalize the skin while reducing the risk of pinpoint bleeding and minimizing potential side effects.

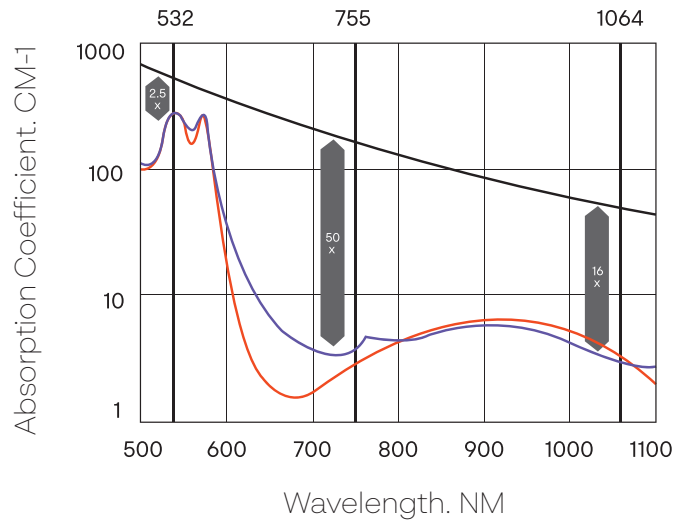

Pre treatment



Immediate post



3 hours post




24 hours post

These before and after photographs are of a patient treated with Cynosure's PicoSure device, not the PicoSure Pro device advertised herein. The PicoSure Pro device has the same parameters as the PicoSure device and so while individual results vary, similar results are expected.

Expert *Voices*/

With more power, new indications, faster treatments and added flexibility, the *PicoSure Pro device takes picosecond technology to unprecedented levels of performance.*

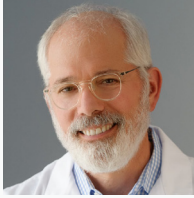

“When it comes to PicoSure Pro, I’ve been really impressed by the power. The additional power allows me to use a bigger spot size and get deeper penetration. I also like the ability to adjust the fluence with different spot sizes. This allows me to customize each treatment to better serve my patients’ specific needs.”

Dr. Emil Tanghetti, Dermatologist



“I love the new PicoSure Pro for treating darker skin types and for treating discoloration. The ability to adjust the fluence for both the flat optic and the platinum focus lens really allows you to finesse each treatment. It’s my device of choice for treating pigment.”

Dr. Nazanin Saedi, Dermatologist



“A man had depression when he visited my clinic for the first time because of the pigmented lesions on his face. After a series of PicoSure treatments, his pigmentation cleared and he is no longer depressed. He is now really positive about his life and I am grateful to have been a part of that.”

Dr. Teppei Sakai, Dermatologist

PicoSure Pro is more than just a laser. It's a tool that empowers you to *transform patients' lives for the better.*

100% patient satisfaction¹

Before & after
3 treatments

Courtesy of Dr. Kentaro Oku



Before & after
3 treatments

Courtesy of Dr. Kentaro Oku



Before & after
3 treatments

Courtesy of Dr. Emil Tangheti



These before and after photographs are of a patient treated with Cynosure's PicoSure device, not the PicoSure Pro device advertised herein. The PicoSure Pro device has the same parameters as the PicoSure device and so while individual results vary, similar results are expected.

Join the *Performance Pro Club*

As a dedicated partner, our goal is to nurture successful partnerships with our customers that go beyond providing best-in-class technology. Through this exclusive club you will gain access to various VIP services, allowing you to focus on what you do best while leaving the rest to us.

*Limited time offer

Exclusive Events

PicoSure Pro events and trainings

Service

Prioritized support

Education & Marketing

Access to the IMCAS Academy E-learning platform as well as a personalized video shoot to utilize on your social platforms*

Networking Events

Access to exclusive VIP lounges and networking events with PicoSure Pro power users

PICOSURE[®] *pro*
by CYNOSURE[®]

